


# How do I select a FermProbe pH electrode?

## FermProbe® pH Electrodes with T-Pull® Cap

F-615 pH electrode with T-Pull cap

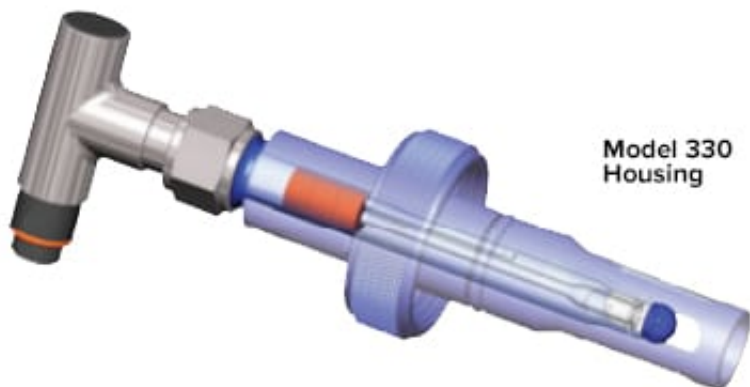

Model F-615 pH FermProbe

The T-Pull cap eases installation and removal of the electrode from the housing and greatly reduces fatigue of the extension cable. The electrode can be removed from the housing without rotating the electrode and twisting the cable, further reducing cable fatigue.

### Features:

- Rugged T-Pull cap eases electrode removal. No tools required.
- T-Pull design greatly reduces cable fatigue.
- Electrode retainer nut is part of handle. Cannot be lost or misplaced.
- Standard S8 detachable cable connection.
- Available in both disconnect cap and integral cabled versions.

---

## FermProbe® pH Electrodes with Standard S8 Metric Cap

### Corporate Office

Broadley-James Corporation  
19 Thomas, Irvine CA, 92618 USA

Phone: 949-829-5555

Tollfree: 800-288-2833

Fax: 949-829-5560

Email: [sales@broadleyjames.com](mailto:sales@broadleyjames.com)

Web: [www.broadleyjames.com](http://www.broadleyjames.com)

### European Office

Broadley-James Ltd.  
Wrest Park, Silsoe  
Beds MK45 4HS, UK

Phone: +44 (0)1525 862518


Fax: +44 (0)1525 862811

Email: [sales@broadleyjames.co.uk](mailto:sales@broadleyjames.co.uk)

Web: [www.broadleyjames.eu](http://www.broadleyjames.eu)

# How do I select a FermProbe pH electrode?

F-635 pH electrode with S8 metric cap

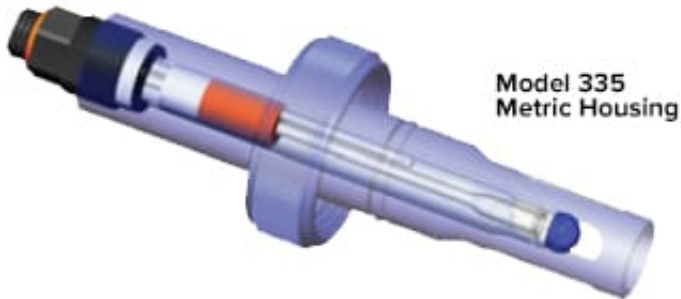

Model F-635 pH FermProbe

This FermProbe style features a standard S8 detachable cable connection on a metric threaded polymer cap that allows the electrode to be used with European style electrode housings. The design requires the rotation of the electrode during installation. It is compatible with Eurodesigned housings requiring a metric Pg13.5 threaded electrode cap.

#### Features:

- Standard S8 detachable cable connection
- Polymer disconnect cap with metric Pg13.5 threads.
- Compatible with European designed bioprocess housings.

---

#### FermProbe® pH Electrodes with K9 Metric Cap

#### Corporate Office

Broadley-James Corporation  
19 Thomas, Irvine CA, 92618 USA

Phone: 949-829-5555  
Tollfree: 800-288-2833  
Fax: 949-829-5560

Email: [sales@broadleyjames.com](mailto:sales@broadleyjames.com)  
Web: [www.broadleyjames.com](http://www.broadleyjames.com)

#### European Office


Broadley-James Ltd.  
Wrest Park, Silsoe  
Beds MK45 4HS, UK

Phone: +44 (0)1525 862518  
Fax: +44 (0)1525 862811

Email: [sales@broadleyjames.co.uk](mailto:sales@broadleyjames.co.uk)  
Web: [www.broadleyjames.eu](http://www.broadleyjames.eu)

# How do I select a FermProbe pH electrode?

F-695 pH electrode with K9 metric cap

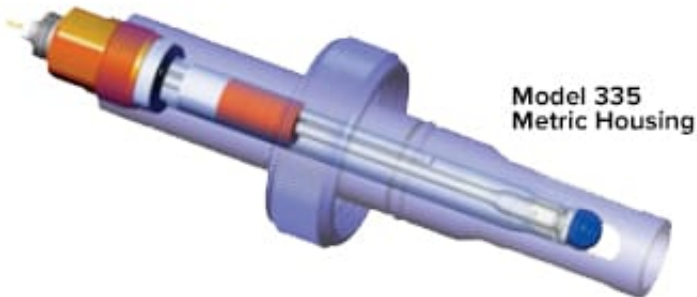

Model F-695 pH FermProbe

This FermProbe style features a K9 detachable cable connection on a metric threaded polymer cap that allows the electrode to be used with European style electrode housings. The design requires the rotation of the electrode during installation. It is compatible with Euro-designed housings requiring a metric Pg13.5 threaded electrode cap.

#### Features:

- K9 detachable cable connection
- Polymer disconnect cap with metric Pg13.5 threads.
- Compatible with European designed bioprocess housings.

#### Corporate Office

Broadley-James Corporation  
19 Thomas, Irvine CA, 92618 USA

Phone: 949-829-5555

Tollfree: 800-288-2833

Fax: 949-829-5560

Email: [sales@broadleyjames.com](mailto:sales@broadleyjames.com)

Web: [www.broadleyjames.com](http://www.broadleyjames.com)

#### European Office

Broadley-James Ltd.  
Wrest Park, Silsoe  
Beds MK45 4HS, UK

Phone: +44 (0)1525 862518

Fax: +44 (0)1525 862811

Email: [sales@broadleyjames.co.uk](mailto:sales@broadleyjames.co.uk)

Web: [www.broadleyjames.eu](http://www.broadleyjames.eu)