


How do I select a FermProbe pH electrode?

FermProbe® pH Electrodes with T-Pull® Cap

The T-Pull cap eases installation and removal of the electrode from the housing and greatly reduces fatigue of the extension cable. The electrode can be removed from the housing without rotating the electrode and twisting the cable, further reducing cable fatigue.

Features:

- Rugged T-Pull cap eases electrode removal. No tools required.
- T-Pull design greatly reduces cable fatigue.
- Electrode retainer nut is part of handle. Cannot be lost or misplaced.
- Standard S8 detachable cable connection.
- Available in both disconnect cap and integral cabled versions.

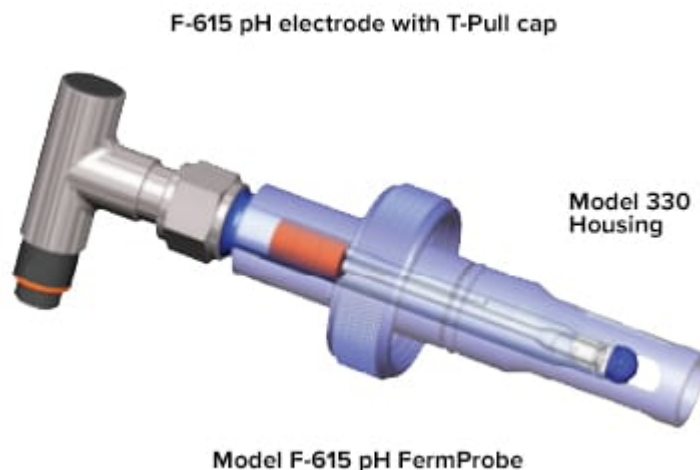

FermProbe® pH Electrodes with Standard S8 Metric Cap

The contents of this publication are presented for information purposes only, and while every effort has been made to ensure their accuracy, they are not to be construed as warranties or guarantees, express or implied, regarding products or services described herein or in their use or applicability. All sales are governed by our terms and conditions, which are available on request. We reserve the right to modify or improve the design or specification of such products at any time.

© 2022 Broadley-James Corporation. All rights reserved. Visit www.broadley-james.com/trademarks for trademark information.

TMP-BF-102101


How do I select a FermProbe pH electrode?

This FermProbe style features a standard S8 detachable cable connection on a metric threaded polymer cap that allows the electrode to be used with European style electrode housings. The design requires the rotation of the electrode during installation. It is compatible with Eurodesigned housings requiring a metric Pg13.5 threaded electrode cap.

Features:

- Standard S8 detachable cable connection
- Polymer disconnect cap with metric Pg13.5 threads.
- Compatible with European designed bioprocess housings.

F-635 pH electrode with S8 metric cap

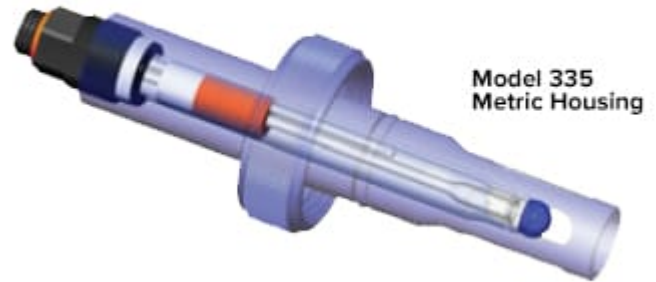

Model F-635 pH FermProbe

FermProbe® pH Electrodes with K9 Metric Cap

The contents of this publication are presented for information purposes only, and while every effort has been made to ensure their accuracy, they are not to be construed as warranties or guarantees, express or implied, regarding products or services described herein or in their use or applicability. All sales are governed by our terms and conditions, which are available on request. We reserve the right to modify or improve the design or specification of such products at any time.

© 2022 Broadley-James Corporation. All rights reserved. Visit www.broadley-james.com/trademarks for trademark information.

TMP-BF-102101


How do I select a FermProbe pH electrode?

This FermProbe style features a K9 detachable cable connection on a metric threaded polymer cap that allows the electrode to be used with European style electrode housings. The design requires the rotation of the electrode during installation. It is compatible with Euro-designed housings requiring a metric Pg13.5 threaded electrode cap.

Features:

- K9 detachable cable connection
- Polymer disconnect cap with metric Pg13.5 threads.
- Compatible with European designed bioprocess housings.

F-695 pH electrode with K9 metric cap

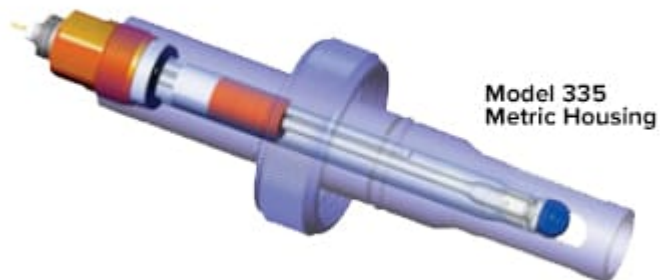

Model F-695 pH FermProbe

The contents of this publication are presented for information purposes only, and while every effort has been made to ensure their accuracy, they are not to be construed as warranties or guarantees, express or implied, regarding products or services described herein or in their use or applicability. All sales are governed by our terms and conditions, which are available on request. We reserve the right to modify or improve the design or specification of such products at any time.

© 2022 Broadley-James Corporation. All rights reserved. Visit www.broadley-james.com/trademarks for trademark information.

TMP-BF-102101